

DrägerService®

**Dräger**

# Field Service Procedure

Part Number: SP00243

Rev: -

Date: 26 April 2002

© 2002 Draeger Medical, Inc.

## **Installation of Auxiliary Fresh Gas Outlet on Fabius GS**

[RETURN TO SERVICE PROCEDURE TABLE OF CONTENTS](#)  
[RETURN TO CD-ROM TABLE OF CONTENTS](#)


## Installation of Auxiliary Fresh Gas Outlet on Fabius GS

### Installation Procedure

**Warning:** Any person involved with the setup, operation, or maintenance of the Fabius GS anesthesia machine must be thoroughly familiar with the instruction manual for this product.


1. Select the correct Fresh Gas label per local application or customer requirement. See Figure 1 for corresponding part numbers.
2. Apply the label to the auxiliary fresh gas outlet where shown in Figure 1. (Both left mount and right mount types are illustrated.)

4117770-001

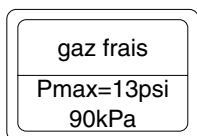


FRENCH

4117770-002

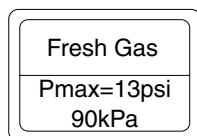


ENGLISH

4117770-004

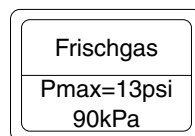


GERMAN

4117770-003

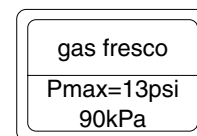


SPANISH

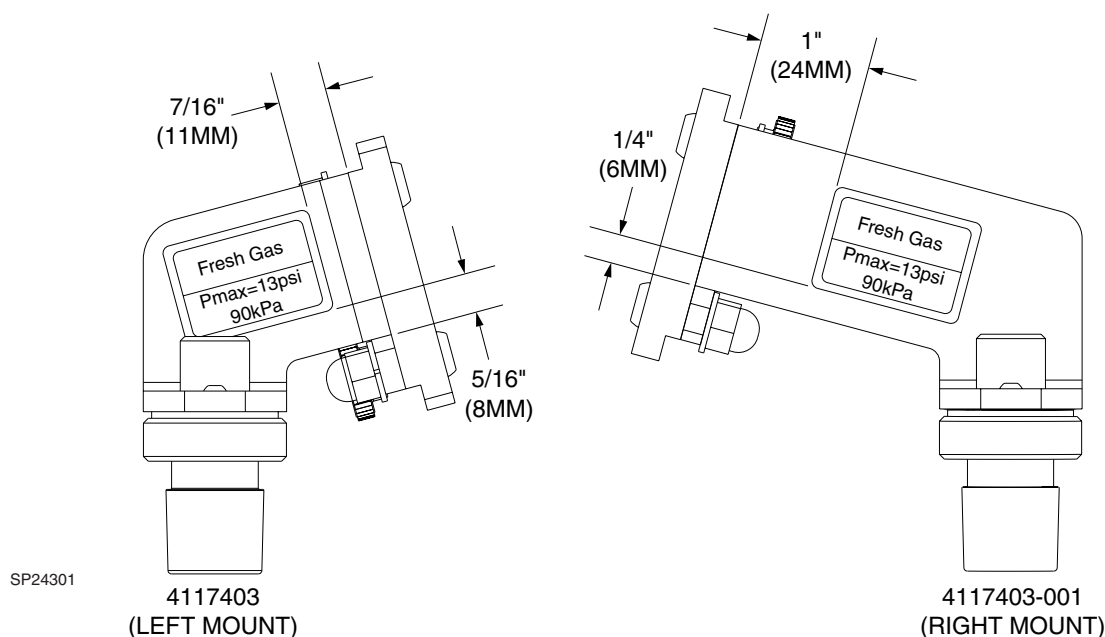

**FIGURE 1. Fresh Gas Label Location  
(approximate measurements)**

## Installation Procedure (continued)


3. Verify the machine is in Standby mode and all flow control valves are closed.
  4. Slide the auxiliary fresh gas outlet assembly up into the appropriate front GCX rail in the machine, and tighten the thumbscrew. (The shear tabs on either side of the slide plate may need to be sanded down to allow the plate to slide freely in the rail.)
  5. Join the machine's fresh gas hose to the rear connector on the auxiliary fresh gas outlet assembly. Left side installation is shown in Figure 2; right side installation is shown in Figure 3.
  6. Connect the Fabius CE hose included in the kit to the auxiliary fresh gas outlet.
- NOTE:** If you are using an A-cone switch, install it on the auxiliary fresh gas outlet and connect the Fabius CE hose to one of the A-cone ports.
7. Connect the other end of the Fabius CE hose to the COSY breathing system.
  8. Pressurize the oxygen system, activate the O2 Flush and verify the fresh gas flow.

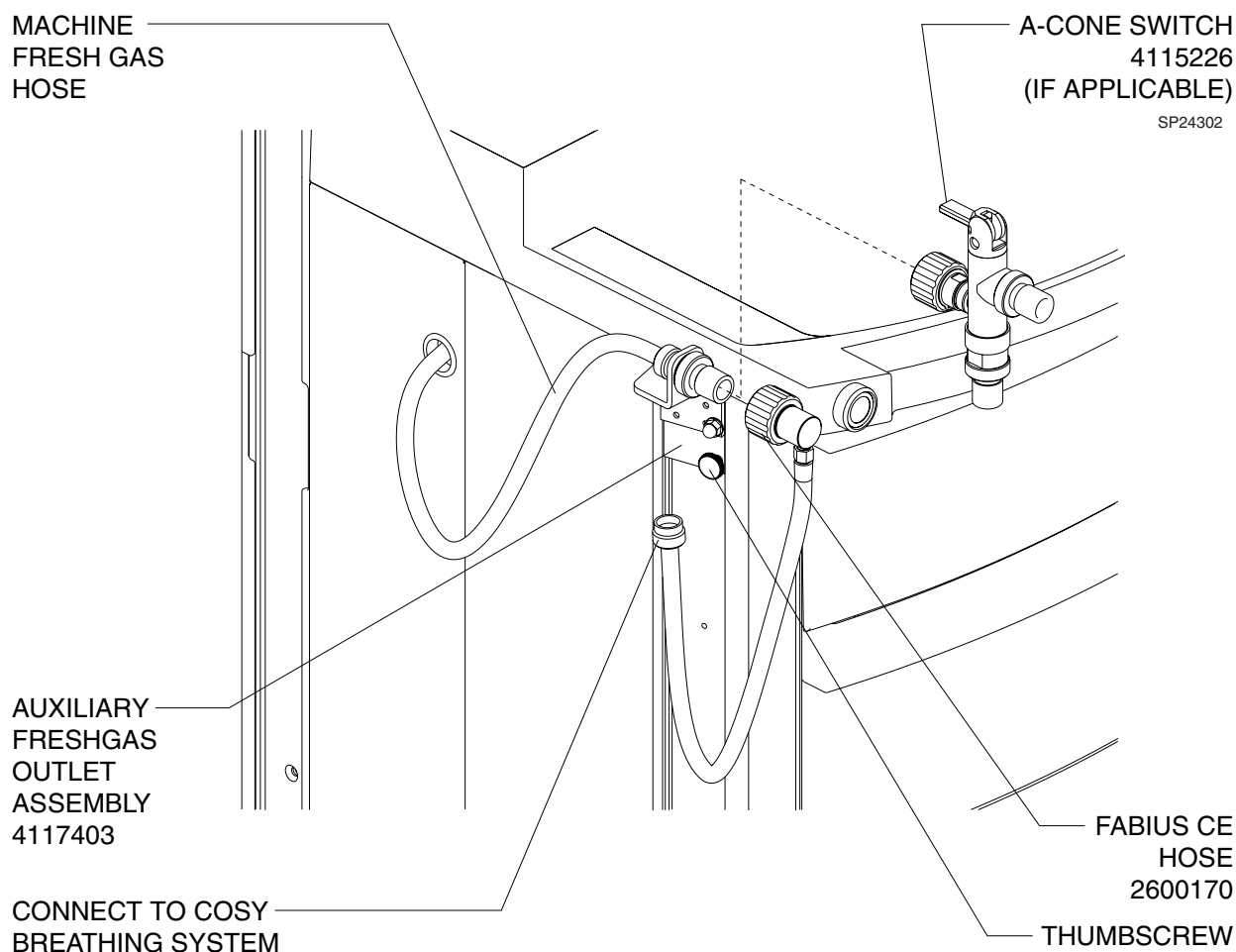

**FIGURE 2. Fresh Gas Outlet Assembly Installation: Left Side**

**Installation Procedure (continued)**

SP24303

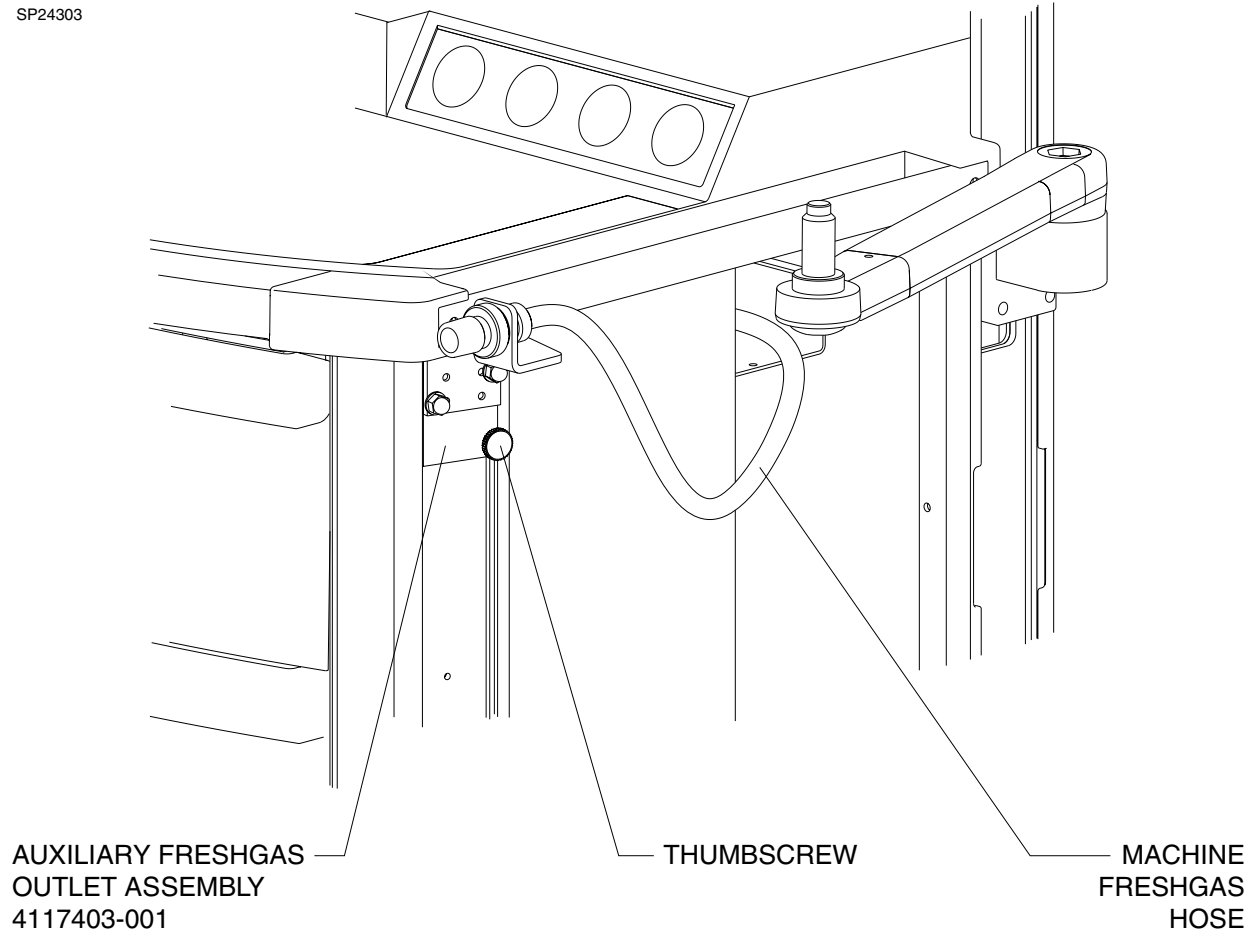

**FIGURE 3. Fresh Gas Outlet Assembly  
Installation: Right Side**

[RETURN TO SERVICE PROCEDURE TABLE OF CONTENTS](#)  
[RETURN TO CD-ROM TABLE OF CONTENTS](#)

[RETURN TO SERVICE PROCEDURE TABLE OF CONTENTS](#)  
[RETURN TO CD-ROM TABLE OF CONTENTS](#)

DrägerService®

**Dräger**

DrägerService is a division of  
Draeger Medical, Inc.  
3122 Commerce Drive  
Telford, PA 18969  
Tel: (215) 721-5402  
(800) 543-5047  
Fax: (215) 721-5784  
Web: [www.draegermedical.com](http://www.draegermedical.com)  
Printed in the U.S.A.